

The r-evolution of the historical model launched in 1978 is now embodied by **Linda Next**.

Linda Next keeps the promise of the fundamental Tecnodent “knee break” concept. The beloved features of compactness, essential design, simplicity and comfort of use are unmistakable.

The dedication to continuous research and development of new technologies and materials brought us to the adoption of RIM (Reaction Injection Molding) system.

The unique result is a combination of flexibility, strength and modern look, with a maximum lifting capacity of 300 Kg.

Extremely intuitive use by one single joystick foot control.

Top performances even in the basic version, including backrest and leg section synchronism, with automatic extensible PU foot section; automatic working programs, table, Trendelenburg, rinse, last and zero positions; soft start and stop and SLOW MODE feature, for precise movements.

Wide, wrapping and slim molded upholstery contributes to enhance the lean image and comfortable feeling.

2 accessory motors for sliding backrest and electrical legrest are just two of the several extras of this chair.

Versatile and comfortable, Linda Next fulfills most demanding needs and requirements.

An outstanding step into the future.

STANDARD AND
OPTIONAL CHARACTERISTICS

• Standard • Optional

	2009next	Lux fit	ECO.next	Linda next	Sting
Headrests	ELLE FISSO	•			
	ELLE/ELLE2 MOBILE	•	•	•	•
Pillows	UNI/Light			•	•
	C95	•	•	•	•
Backrest	C2002	•	•	•	•
	Backrest motorized height adjustment				
	Narrow		•	•	
	Large	•	•	•	•
Armrest	“T” shaped				
	Right armrest	•	•	•	•
	Left armrest	•	•	•	•
Programs	3 customizable programs (T2000 PCB)	•	•	•	•
	Automatic built in programs: Patient Positioning, Table, Trendelenburg				•
	9 user programmable positions (3x profile)				•
	Return to “zero” program	•	•	•	•
	Rinse program	•	•	•	•
Unit connection	Double side drilled upper body	•	•	•	•
	Fixed central arm	•	•	•	•
	Rotation	•	•	•	•
	Rotation with central arm	•	•	•	•
Function	Rear connection for unit bracket	•	•	•	•
	Trendelenburg				•
	Knee-break legrest				•
	Independent legrest movement				•
	Securities	•	•	•	•
Control	Small integrated utility box	•			
	Large integrated utility box	•	•		
	Joystick	•	•	•	•
	Mobile foot control	•	•	•	•
	Infrared remote control	•	•	•	•

HEADRESTS

Tecnodent manual Headrests are available in 4 different models: ELLE FISSO, ELLE/ELLE2 MOBILE, UNI and Light. Every headrest was made for meeting every operator’s need.

ELLE/ELLE2. ELLE/ELLE2 FISSO has only the movement along its vertical axis. ELLE/ELLE2 MOBILE is an anatomical headrest, characterized by movement of inclination of the upper part (double articulation) for a correct cervical support.



ELLE/ELLE2



C95



C2002

Magnetic pillows. ELLE/ELLE2 type, both fix and mobile version, can be combined with magnetic pillow C95 or C2002. Either pillow can be applied directly on the backrest, for ensure a better comfort, independently of the patient’s height.



UNI



Light

UNI and Light. The UNI and Light headrests are multi-articulated and allows positioning of the headrest along 3 axis.



Tecnodent Exim India (P) Ltd.

B138, Ground Floor, Mohammadpur, Bhikaji Cama Place
New Delhi , India - 110066

Phone: +91 11 41550193, +91 99102 11888 | Email: TEI@TEIDental.com

Web: www.teidental.com



LINDA NEXT STANDARD FEATURES

- Complete symmetry for right – or left-handed professionals;
- Electromechanical lifting, lowering and backrest tilting;
- Backrest movement synchronized with leg section;
- Automatic extensible PU foot section;
- Independent Trendelenburg movement;
- Automatic built in programs for patient positioning: table, Trendelenburg, rinse, last and zero position;
- 3 operator profiles with 3 programmable positions each;
- Safety security STOPS on backrest, legrest and controls;
- SLOW MODE feature for precise movements;
- Maximum lifting capacity 300 Kg;
- 24 DC linear actuators;
- Soft start & soft stop;
- Joystick control;
- ELLE FISSO headrest;
• Standard molded upholstery;
- RAL 9003 painting;
- Different PU and fabric colors.



TECHNICAL DATA

Minimum height	mm	500
Maximum height	mm	875
Backrest height adjustment	mm	60
Backrest tilting angle	°	0°-90°
Width	mm	630
Length	mm	2000
Net weight approx.	Kg	140
Maximum lifting capability (patient + unit)	Kg	300
Lifting raising speed	mm/sec	22
Trendelenburg	°	21,5°
Rotation angle	°	325°
Classification		I
Class		B

LINDA NEXT - Options


- Backrest motorized height adjustment
- Backup Battery
- Children pillow
- Different chair painting colors
- Different headrests available according to specific treatment needs
- European or American unit bracket
- Independent leg section control
- IR remote control
- Leg section foot protection
- Memory or hand made soft foam upholstery
- Rotation
- Stabilizing plate
- Swiveling Left & Right armrests
- Wheels




Motorized backrest height adjustment



Swiveling armrests

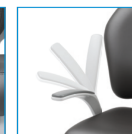

IR remote control

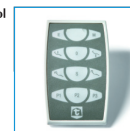

Soft Foam

