

The latest evolution of the unmistakable and historical Linda 3 model has arrived: **LINDA EVO** is finally here.

**LINDA EVO** is the ultimate mid-range knee-break chair, perfectly suitable to treat all patients, even children, elderly people, and persons with reduced mobility, for both dental and medical procedures.

The distinctive features of the well renown Tecnodent knee-break chair concept, such as compactness, easy accessibility, clean lines, strong solidity, and comfort of use, are preserved and empowered with modern design, technology, and structure.

Low voltage motors, with Hall sensor limit switches, guarantee smoothness and perfect movement fluidity, thanks also to soft-start and soft-stop features, while the unmistakable designer armrests and split-section upholstery, along with the enveloping backrest, uphold the smart design statement.

Versatile and comfortable, **LINDA EVO** embodies the perfect “multifunctional smart knee-break”, thanks to its elegant yet compact, ergonomic, and functional style.

STANDARD AND OPTIONAL CHARACTERISTICS		LINDA NEXT	LINDA EVO
	• Standard • Optional		
Headrests	ELLE FISSO	•	•
	ELLE/ELLE2 MOBILE	•	•
	UNI/Light	•	•
	Vale		•
Pillows	C95	•	•
	C2002	•	•
Backrest	Backrest motorized height adjustment	•	
	Large enveloping	•	
	High enveloping		•
Armrest	Right armrest	•	•
	Left armrest	•	•
Programs	Automatic built in programs: Patient Positioning, Table, Trendelenburg	•	
	9 user programmable positions (3x profile)	•	
	3 user programmable positions		•
	Return to “zero” program	•	•
	Rinse program	•	•
Other features	Rotation	•	•
	Lamp / unit bracket	•	•
	Back-up Battery	•	
	Wheels	•	
Function	Trendelenburg	•	•
	Knee-break legrest	•	•
	Independent legrest movement	•	
Control	Securities	•	•
	Joystick	•	•
	Infrared remote control	•	•

**HEADRESTS**

Tecnodent manual Headrests are available in 5 different models: ELLE FISSO, ELLE/ELLE2 MOBILE, UNI, Light and VALE. Each headrest was made for meeting every operator’s need.

**ELLE/ELLE2.** ELLE/ELLE2 type, both fix and mobile version, can be combined with magnetic pillow C95 or C2002. Either pillow can be applied directly on the backrest (standard or memory-foam only) for ensuring a better comfort, independently of the patient’s height.



ELLE/ELLE2



C95



C2002

**Magnetic pillows.** ELLE/ELLE2 type, both fix and mobile version, can be combined with magnetic pillow C95 or C2002. Either pillow can be applied directly on the backrest, for ensure a better comfort, independently of the patient’s height.



UNI



Light

**UNI and Light.** The UNI and Light headrests are multi-articulated and allow positioning of the headrest along 3 axis.



VALE

VALE is a double articulated headrest with enveloping cushion.





## LINDA EVO STANDARD FEATURES

- Fully electromechanical upward, downward and tilt movements of the backrest
- Backrest movement compensated, up to horizontal position
- Tilting of the backrest synchronized with the elongation of the leg rest
- Headrest model ELLE FISSO (vertical adjustment only)
- High enveloping backrest
- Complete symmetry, usable by both right and left-handed professionals
- Safety stops on backrest and lifting arm
- Automatic program for complete return to the patient lowering position ("zero")
- Three user programmable customized working positions
- Automatic rinse program with return to the last working position
- Low voltage controls: 5 V
- Joystick control, on the base
- Soft-start and soft-stop movement



### TECHNICAL DATA


Max. height	mm	880
Min. height	mm	550
Backrest tilting angle	°	0 -76
Seat tilting angle	°	+7/+18
Width	mm	615
Total length	mm	1750
Gross weight	Kg	110
Net weight	Kg	80
Maximum lifting capability (patient + connections)	Kg	135+75
Trendelenburg	°	+18
Rotation angle	°	+/-160
Classification		I
Cass		B

## LINDA EVO - Options

- Double-articulated headrest model VALE
- Double-articulated headrest model ELLE MOBILE
- Head rest model ELLE-2 MOBILE
- Magnetic cushion C95
- Magnetic cushion C2002
- Multi-articulated headrest model "UNI tappezzato"
- Multi-articulated headrest model "LIGHT"
- Right and/or left armrest
- Electronic board with lamp activation relay
- Rotation system around the vertical axis under the chair base
- Infrared hand remote control
- Children booster seat
- ENT footrest
- Seamless memory foam upholstery
- Handmade soft-foam upholstery
- Handmade soft-foam upholstery with contrasting color stitching




Joystick pedal control

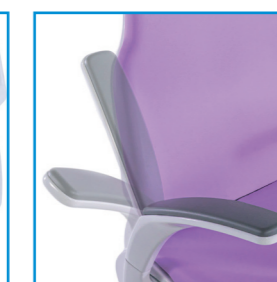

Swiveling armrests

IR remote control

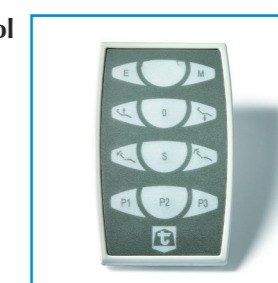

Soft Foam